

# MEDIOSTAR

## NEXT FAMILY



THE UNMATCHED STARS FOR  
**PROFESSIONAL HAIR REMOVAL**

Brochure not for the U.S.

  
**Asclepion**  
Laser Technologies



THE NEW  
**MEDIOSTAR NeXT FAMILY**  
BEYOND PROFESSIONAL EPILATION

**MEDIOSTAR**  
**NeXT XL**

The largest spot size available on the market and treatments faster than ever.  
Try today the power of an extraordinary hair removal experience. >>>

**MEDIOSTAR**  
**NeXT ALEX**

The efficiency of the Alex wavelength meets the latest diode technology.  
A unique combination that allows professional results even on thinner  
and lighter hairs. >>>

**MEDIOSTAR**  
**NeXT PRO**

The best choice for users who are starting in the hair removal business  
and want to get on board with the latest technology.  
Reliable, with no running costs, upgradable: simply ready for the success. >>>

**MEDIOSTAR**  
**NEXT XL**



**XTRAPOWERFUL**  
**XTRAFAST**  
**XTRALARGE**

# FASTER THAN EVER

The new MeDioStar NeXT XL features the largest spot available on the market and offers an extraordinary hair removal experience at the highest speed. Thanks to **SmoothPulse Mode** the treatment is faster than ever and still painless. The new **360° skin cooling system** allows a highly efficient pre-cooling of the epidermis, protecting it from the risk of burns and allowing the treatment in any direction. This guarantees the best comfort for both the patient and the user. The **XL-S and XL-L tips** are specially designed for the treatment of both larger and smaller areas.

All the flexibility that you need in an XL size.

+ LARGEST SPOT ON THE MARKET

+ TREAT ALL SKIN TYPES,  
ALL YEAR ROUND

+ SMOOTHPULSE MODE FOR  
PAINLESS TREATMENTS

+ 360° SKIN COOLING  
SYSTEM INTEGRATED

+ NO CONSUMABLES

Spot size: XL-L 9.1 cm<sup>2</sup>

Spot size: XL-S 3.0 cm<sup>2</sup>



TREAT A BACK IN  
LESS THAN 4 MINUTES



# MEDIOSTAR

NEXT ALEX



**ALEX WAS GOOD.**  
**WE MADE IT BETTER.**

# THE NEXT ALEX GENERATION

The MeDioStar NeXT ALeX brings the **Alexandrite wavelength and the latest Diode Technology** together, combining the efficacy of the high melanin absorption of the Alex with all the benefits of the Diode. Thanks to the highly selective 755 nm wavelength, even thinner and lighter hairs can be treated with the best results. The **SmoothPulse Mode** allows fast and painless treatments and the **integrated skin cooling system** guarantees the highest safety. The proven reliability of the Diode Technology assures **the optimization of the treatment's procedure** (no time consuming heating up) and **no running costs**.

Recovering your investment has never been so easy.

+ HIGHLY EFFECTIVE ON THINNER AND LIGHTER HAIRS

+ TREAT UP TO SKIN TYPE IV

+ INTEGRATED SKIN COOLING SYSTEM

+ SMOOTHPULSE MODE FOR PAINLESS TREATMENTS

+ FASTEST RETURN ON INVESTMENT

Spot size: 1.4 cm<sup>2</sup>



**FINALLY THE ALEX  
WAVELENGTH WITH ALL  
THE BENEFITS OF THE DIODE**

**MEDIOSTAR**  
NEXT PRO

**THE UPGRADABLE SOLUTION.**  
PREPARE FOR SUCCESS.





## GET ON BOARD WITH THE LATEST TECHNOLOGY

The MeDioStar NeXT PRO is the most suitable solution for new users who want to start the hair removal business with the state-of-the-art technology. Thanks to the **highly reliable Diode Technology** and the **user-friendly interface**, even less experienced users can become confident with the most advanced hair removal techniques in a few simple steps. The **SmoothPulse Mode** and the **integrated skin cooling system** allow fast, safe and effective treatments, with basically no pain. On top of this, the device allows **multiple treatments** including **acne treatment**, **skin rejuvenation** and **skin whitening** and it is easily upgradable also for **vascular treatments**.

It is just exactly what you need for getting off to a good start.

+ FAST, SAFE AND EFFECTIVE

+ TREAT ALL SKIN TYPES,  
ALL YEAR ROUND

+ INTEGRATED SKIN COOLING SYSTEM

+ SMOOTHPULSE MODE FOR  
PAINLESS TREATMENTS

+ MULTIPLE TREATMENT OPTIONS

Spot size: 1.4 cm<sup>2</sup>



**UPGRADABLE & FLEXIBLE,  
THE BEST SOLUTION FOR  
YOUR BUSINESS**

# UPGRADES\*

## MAKE THE MOST OF YOUR DEVICE

---

### BEYOND JUST HAIR REMOVAL

## VASCULAR HANDPIECE

MeDioStarNeXT goes beyond just hair removal, with versatility and user-friendliness in other applications, such as the **treatment of leg veins**, hemangiomas and other ectatic vessels. With the special vascular handpiece featuring a wavelength of 940 nm, it allows non-invasive and fast treatments, usually with no need for an anesthetic.



**Spot size: 0.12 cm<sup>2</sup>**

### POWER IN YOUR HANDS

## HIGH POWER HANDPIECE

Patients with **darker skin types and very fine hairs** can be very difficult cases to deal with. The High Power handpiece offers the possibility to treat this condition with highest effectiveness and completely safely, thanks to the diode technology reaching up to 90J/cm<sup>2</sup> with very short pulses.



**Spot size: 0.48 cm<sup>2</sup>**

### THE PORTABLE SOLUTION

## TROLLEY

The specially designed trolley allows the easiest transportation and the **maximum space optimization**, including a storage space for footswitch and power cord in the back. A stylish additional feature, for those who like the maximum comfort.




\*optional

# ASCLEPION

## COMPETENCE, EXPERIENCE, SUCCESS

Asclepion Laser Technologies has been operating for **more than 35 years** as a leader on the international medical laser scene as manufacturer of advanced laser systems for dermatology, aesthetic medicine and surgery. Based in the Jena Optical Valley, famous throughout the world as the homeland of the German optics industry, the company invests continuously in research and development of new technologies. Today, customers in **more than 60 countries** trust in the Made in Germany technology of Asclepion and its proven scientific expertise.

**MeDioStar NeXT, the 5<sup>th</sup> generation of Asclepion's hair removal technology, boasts an unmatched success all around the world, proof of its effectiveness, safety and high quality.**

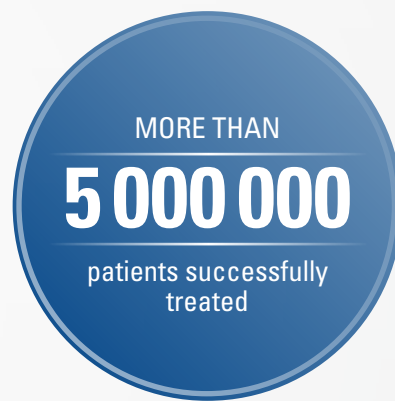

High quality means for us also a customer oriented philosophy. That's why we work hard every day in order to provide you not only the best technology but also a full range of services, to support you at 360°.

### TRAINING AND EDUCATION



Both regional and international trainings are permanently organized in our training center in Jena. That's how we provide our customers all the knowledge that they need for safely and effectively using our laser devices.

Find out more at:  
[www.asclepion.com/academy](http://www.asclepion.com/academy)

### MARKETING AND WEBCLUB



A wide range of marketing tools is available on the Asclepion's WEBCLUB, an online platform designed for giving you a real time access to all the latest news and important documents for the success of your practice.

Find out more at:  
[www.asclepion.com/webclub](http://www.asclepion.com/webclub)

### TECHNICAL SERVICE



Our highly specialized technical assistance, provided both directly and through a network of local distributors, follows the customer from the installation, to the assistance in case of defect, to the delivery of spare parts.

Find out more at:  
[www.asclepion.com/service](http://www.asclepion.com/service)



# Asclepion

## Laser Technologies

### SPECIFICATIONS

<b>Laser</b>	High Power Diode Laser, class 4
<b>Power</b>	Up to 2400 W
<b>Wavelength</b>	755–950 nm*
<b>Fluence</b>	Up to 90 J/cm <sup>2</sup> (hair removal)* Up to 210 J/cm <sup>2</sup> (vascular treatments)
<b>Spot size</b>	9.1, 3.0, 1.4 or 0.48 cm <sup>2</sup> (hair removal)* 0.12 cm <sup>2</sup> (vascular treatments)**
<b>Frequency</b>	Up to 12 Hz*
<b>Display</b>	8.4" LCD Display
<b>Sizes</b>	Table Top: 36.5 x 56 x 46.5 cm With Trolley**: 36.5 x 56 x 93 cm
<b>Weight</b>	Approx. 35 kg

\* according to the handpiece    \*\* optional  
All specifications are subject to change without notice



Patents  
DE 197 10 676  
US 6,162,218  
Not for sale in the US  
Certified according to  
DIN EN ISO 9001: 2008  
DIN EN ISO 13485: 2010  
Medical Device Directive  
93/42/EEC, Annex II

 **MADE IN GERMANY**

FOLLOW US ON

